


SealBio-1

Automated Sealer

SealBio-1 automated sealer is a device for rapid and continuous sealing of microplates. It can be adapted to ANSI, SBS standard 12-well, 24-well, 48-well, 96-well and 384-well plates with height range of 3~62 mm, including all kinds of PCR plates, ELISA plates, deep-well plates, cell culture plates, etc. It supports the use of 78 mm puncturable roll film, aluminum foil roll film and strong heat sealing roll film.

Suitable for pre-filled of biochemical reagents, sealing film before PCR assays, long-term storage of samples, and other experimental applications that require sealing of reagents.

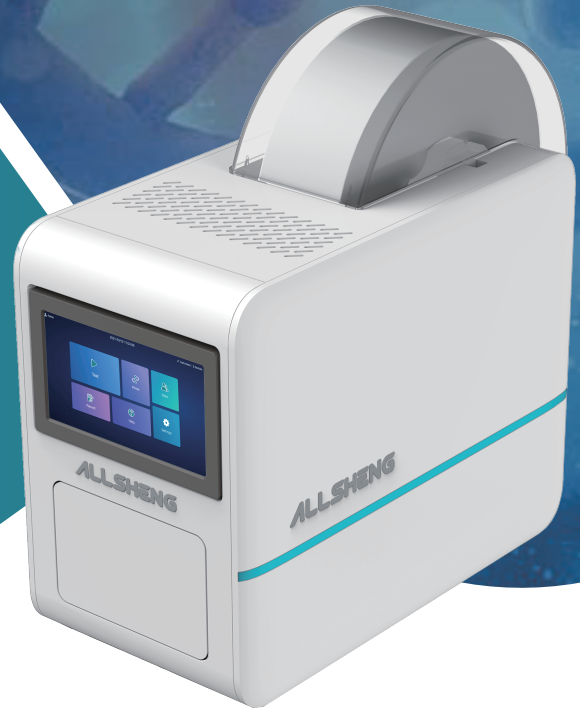

✓ RT.~200 °C free setting, to meet the needs of different film sealing applications

✓ Built-in independent air source, no need to connect other equipment

✓ Automatic preloading of film, continuous operation

Film sealing cycle time <20 seconds ✓

Touch screen operation, graphical interface is easier to use ✓

Standard roll film protection cover to avoid film surface contamination ✓

Program storage function to save different types of film parameters for quick recall ✓

Specification

Parameter / Model	SealBio-1	
Product name	Automated sealer	
Supported plates	ANSI / SBS standard plate, height range 3~62 mm	
Support roll film	Film type	Puncturable film / aluminum foil film / strong heat sealing film
	Film width	78 mm
	Loading method	Automatic preloading of roll film
Film sealing temperature	Temperature range	RT.~200 °C
	Temperature control accuracy	±1 °C
	Accommodative amplitude	1 °C
Pre-heating time	5 min	
Cycle time	<20 seconds/time	
Adjustable film sealing time	Time range	0.5~9.9 seconds
	Accommodative amplitude	0.1 seconds
Air source method	Built-in air source	
Port	1×USB-B port and 1×Ethernet port	
Power supply	AC 100~240 V, 50/60 Hz	
Maximum power	300 W	
Software operation	7-inch touch screen software (resolution 1024*600) / edit sealing scheme / prestored sealing scheme	
Safety	Standby function (low power mode) / detection of plate /detection of roll film	

Ordering Information

Code	Description
AS-15010-00	SealBio-1 automated sealer

HANGZHOU ALLSHENG INSTRUMENTS CO., LTD.

Building 9 No.7 of Zhuantang Science and Technology Economic Zone,
Xihu District, Hangzhou City, 310024 Zhejiang, P.R. China

Tel: +86-571-88859758

Fax: 0086-571-87205673

✉ info@allsheng.com

🌐 www.allsheng.com

