

Smart Frequency Conversion ULT Freezer



Core Advantages

- IoT technology for more intelligent management
- HC & Inverter technology for superior energy savings

Friendly Design

- Safe and secure
- Cloud data storage available
- Low noise design, reducing the noise down to 53dB
- Optimized insulation

Qingdao Haier Biomedical Co.,Ltd.

No.280 Feng Yuan Road, High-tech Zone,
Qingdao, 266109, P.R. China
Tel: +86-0532-88935593
Website: www.haiermedical.com



Haier Biomedical
International



Haier Biomedical
International



@haiermedicalint



Haier Biomedical
International



Haier Biomedical
International

Scope of Application



Applicable for products and samples which require strict storage conditions such as viruses, pathogens, red blood cells, white blood cells, skin, bones, bacteria, semen, biological products, electronics and special materials. Suitable for long term storage applications and compliant with typical storage requirements found in hospitals, disease control and prevention centres, scientific research institutions biomedical engineering institutes, agriculture/ fishery companies as well as the electronics and chemical industry.



DW-86L829BPT


Smart Frequency Conversion ULT Freezer

Advanced Hardware System



Smart Full-size Touch Screen

10 inch capacitive touch screen, state-of-art UI design coupled with sample management system to provide the optimal user experience.

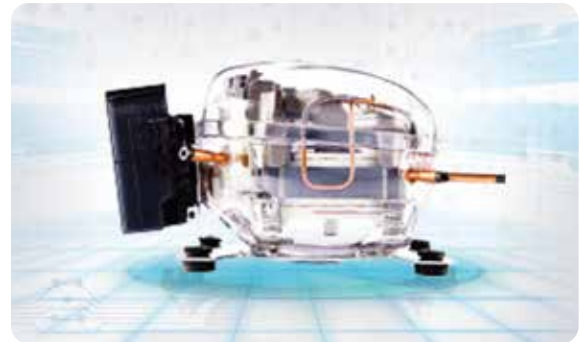

HC & Variable Frequency Drive Refrigeration System for Additional Energy Saving

Advanced innovative design delivers superior energy savings. The energy consumption is down to a single digit.

Optional IoT Software System



Simplified Sample Management Experience

Optional barcode scanner for simple, effortless and precise identification. Input and retrieve your samples with higher precision and efficiency.



Wireless Monitoring Connectivity

Check the real-time operating status via mobile phones or portable tablet, simple and reliable.

A Quicker product access, identification and retrieval

Instead of a manual system, the one-gun, one-code and one-key operation plus touch screen synchronization means you can access and retrieve your samples within seconds.

B 24-hour sample protection

Using Haier Biomedical's app and IoT technology the unit can be monitored and can self-diagnose faults, ensuring you are always aware of your unit's status and able to make real-time and informed choices to protect your samples.

C Cloud data storage available

Store hundreds of millions of scientific research and sample information in the cloud server.

Friendly Design



Safe and secure

Standard equipped with key lock, padlock and electromagnetic lock with optional fingerprint lock, providing multiple safeguards for sample safety.



Low noise design, reducing the noise down to 53dB

Special noise-reduction design plus super silent compressor technology and energy-saving fan, considerably lowers noise level.



Super energy efficient with three environmentally friendly innovations

Ultra-low power consumption, down to less than 10kWh/day, ensuring a world leading energy saving performance.



Optimized insulation

Double foaming for both inner and outer doors and five-layer sealing design and optimized super-thick VIP thermal insulation technology, extends temperature holdover time during power failure and increases insulation efficiency by 20%.

Specifications



Model		DW-86L579BPT	DW-86L729BPT	DW-86L829BPT	DW-86L959BPT	
Technical Data	Cabinet Type	Upright	Upright	Upright	Upright	
	Climate Class	N	N	N	N	
	Cooling Type	Direct Cooling	Direct Cooling	Direct Cooling	Direct Cooling	
	Defrost Mode	Manual	Manual	Manual	Manual	
	Refrigerant	HC	HC	HC	HC	
	Sound Level (dB(A))	43.5	43.5	43.5	47	
Performance	Cooling Performance (°C)	-86	-86	-86	-86	
	Temperature Range (°C)	-40~-86	-40~-86	-40~-86	-40~-86	
Control	Controller	Microprocessor	Microprocessor	Microprocessor	Microprocessor	
	Display	LCD Touch Screen	LCD Touch Screen	LCD Touch Screen	LCD Touch Screen	
Electrical Data	Power Supply (V/Hz)	100~230/50/60	100~230/50/60	208~230/50/60	208~230/50/60	
	Power (W)	1100	1100	1100	1300	
	Electrical Current (A)	6	6	6	7	
Construction	Capacity (L/Cu.Ft)	579/20.4	729/25.7	829/29.2	959/33.9	
	Net/Gross Weight (approx)	kg	325/355	350/385	380/415	450/485
		lbs	716.5/782.6	771.6/848.8	837.7/914.9	992.1/1069.2
	Interior Dimension (W*D*H)	mm	620*716*1310	766*716*1310	870*716*1310	1016*716*1310
		in	24.4*28.2*51.6	30.2*28.2*51.6	34.3*28.2*51.6	40.0*28.2*51.6
	Exterior Dimension (W*D*H)	mm	895*998*1980	1046*998*1980	1145*998*1980	1296*998*1980
		in	35.2*39.3*78.0	41.2*39.3*78.0	45.1*39.3*78.0	51.0*39.3*78.0
	Packing Dimension (W*D*H)	mm	950*1055*2150	1100*1105*2150	1190*1045*2150	1365*1105*2150
in		37.4*41.5*84.6	43.3*43.5*84.6	46.9*41.1*84.6	53.7*43.5*84.6	
Loading Quantities	Container load (20'/40'/40'H)	12/24/24	10/20/20	8/20/20	8/16/16	
Alarms	High/Low Temperature	Y	Y	Y	Y	
	Hot Condenser	Y	Y	Y	Y	
	Power Failure	Y	Y	Y	Y	
	High/Low Voltage	Y	Y	Y	Y	
	Sensor Error	Y	Y	Y	Y	
	Low Battery	Y	Y	Y	Y	
	High Ambient Temperature	Y	Y	Y	Y	
	Door Ajar	Y	Y	Y	Y	
Accessories	Caster	Y	Y	Y	Y	
	Foot	Y	Y	Y	Y	
	Porthole	Y/2	Y/2	Y/2	Y/2	
	Shelves/Inner Doors	3/4	3/4	3/4	3/4	
	USB Interface	Y	Y	Y	Y	
	Remote Alarm (Dry contacts)	Y	Y	Y	Y	
	5V Power Supply Port	Y	Y	Y	Y	
	Temperature Recorder	Optional	Optional	Optional	Optional	
	RS485 Port	Y	Y	Y	Y	
	CO ₂ Backup System	Optional	Optional	Optional	Optional	
	LN ₂ Backup System	Optional	Optional	Optional	Optional	
Certifications	CE	Y	Y	Y	Y	
	UL	Y	Y	Y	Y	
	ENERGY STAR	Y	Y	Y	Y	

Product appearance and specifications are subject to change without notice