

Pharmacy Refrigeration

Portable Cooler Range

Portable Vaccine Cooler 18L

GDPA code **VCP18UK 18L Portable Cooler**
 SKU code **444410883**




Key features

- ✓ +2°C to +8°C operating temperature range
- ✓ External LED temperature display
- ✓ Holds temperature inside between 2°C & 8°C for up to 30 hours
- ✓ Aluminium heat preserving foil lining
- ✓ Typical vaccine storage 12-13L
- ✓ Covered bung non slip feet
- ✓ Webbing adjustable shoulder strap
- ✓ Comes with 12 gel packs

Specification: Portable Vaccine Cooler 18L

Height (mm)	310
Width (mm)	440
Depth (mm)	285
Weight (kg)	4.03
Gross Capacity (L)	18
Net Capacity (L)	8
Temperature Range	2 to 8 degrees celsius
Ambient Temperature Range	-10 to +40 degrees celsius
Temperature Hold Range (hours)	Upto 30
Aluminium Lined Sleeve	Yes
Handle Type	Plastic
Additional Shoulder Strap	Yes
Gel Packs Included	Yes
Number of Gel Packs	12
CE Certified	Yes
External Temperature Display	Yes
Internal Air Temperature Probe	Yes
Battery Included	Yes
Lockable	No

18L Manual Portable Cooler - 25C / 50% rh - 290719

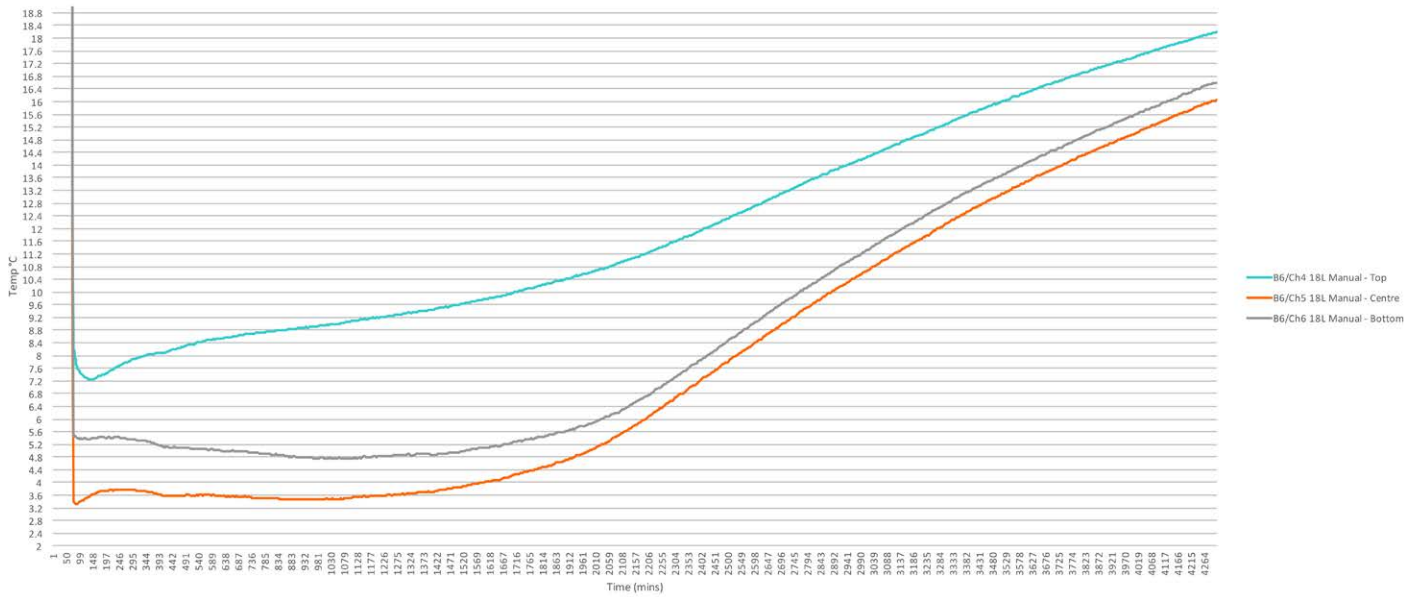

18L manual portable cooler:

Time to reach above 8°C = 42 hrs 14 mins
1.76 Days
(2536 mins)

Load Set Up

Ice packs in freezer at -18.9-19.6°C for 64 hrs
Gel packs in Ward fridge (1-5C) - display settles at 2.3°C for 42.5 hrs
Fully loaded with Gel packs