tel. (+48) 32 453 91 70

A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

Data sheet

Cooled Incubator (IL) **ILW 750 Smart PRO**



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART PRO controller: - large (7"), clear, full colour touch screen

- LAN, USB ports and WiFi for communication and data transfer
- multi-segment time and temperature programs
- overview of data in tabular and graphic form
- visual and sound alarm
- Admin function for management
- password protected log-in
- internal memory for programs and data storage
- event registry with user notifications
- LabDesk software and user manual for direct download



Smart PRO - preview screen






A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

| TECHNICAL DATA | |
|--|---|
| air convection | forced |
| chamber capacity [I] | 749 |
| working capacity [I] | 749 |
| controller | microprocessor PID |
| display | 7" full colour touch screen |
| TEMPERATURE | |
| temperature range [°C] | 0+100 / -10+100 (option) |
| temperature resolution every [°C] | 0,1 |
| temperature fluctuation at 37°C [+/-°C]* | 0,2 |
| temperature variation at 37°C [+/-°C]* | 0,3 |
| over temperature protection | class 3.3 to DIN 12880 |
| CHAMBER | |
| door type | double(5) / door with viewing window (option) |
| interior | |
| Smart PRO | DIN 1.4301 |
| IG Smart PRO | DIN 1.4301 |
| housing | |
| Smart PRO | powder coated sheet |
| IG Smart PRO | stainless steel linen finish |

| overall dims [mm] /1/ | |
|-----------------------|------|
| width A | 1260 |
| height B | 1910 |
| depth C | 890 |
| internal dims [mm] | |
| width D | 1040 |
| height E | 1200 |
| depth F | 600 |

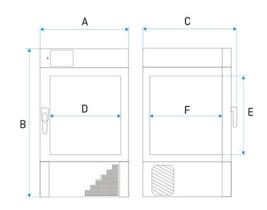

| shelves (standard max) | 5 16 |
|--|--------|
| - reinforced shelf version (PW) [kg] /3/ | 100 |
| max unit workload [kg] | 140 |
| - reinforced unit version (W) [kg] /4/ | 300 |
| weight [kg] | 256 |





A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

| ELECTRICAL PARAMETERS | | |
|-----------------------|-------------------|--|
| voltage** | 230V 50-60Hz | |
| nominal power [W] | 1900 | |
| refrigerant | R290 / GWP=3 | |
| warranty | 24 months | |
| manufacturer | POL-EKO-APARATURA | |

all the above technical data refer to standard units (without optional accessories)

st - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

K = +/- (T average max. - T average min.) / 2

- ** other power supplies on request
- 1 depth doesn't include 50 mm of power cable
- 2 on uniformly loaded surface
- 3 reinforced shelf
- 4 additional internal glass door

OPTIONS AND ACCESSORIES

| | Order number: */A | door with viewing window |
|------|---------------------------|---------------------------------|
| INOX | Order number: */P INOX | Stainless steel wire shelf INOX |
| | Order number: */PP | Perforated shelf |
| | Order number: */PW | Reinforced shelf |
| | Order number: KUW GN*/* | Stainless steel cuvettes |
| *** | Order number: *PLUS | Automatic defrosting function |
| * | Order number: */T | Low temperature version |
| KG | Order number: */W | Reinforced version |
| | Order number: OWW/OWW LED | Interior lighting |
| | Order number: */FIT | Phytotron system |



tel. (+48) 32 453 91 70



A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

| 90% 30% | Order number: FIT/R3 | FIT panels independent control |
|------------|-----------------------------------|--|
| • | Order number: PT 100 | Additional temperature sensor |
| | Order number: GNZ | Internal socket |
| >0< | Order number: BRT/*/L or IQ/OQ/PQ | Calibration and IQ, OQ, PQ qualification |
| 4 | Order number: BPP 12 | Battery backup for display |
| ALARM | Order number: PORT ALARM | Dry alarm contact |
| USB | Order number: USBK | USB cable |
| LAN | Order number: LANK | LAN cable |
| 20mm | Order number: OCZ/20 | Non-standard access port 20 mm |
| 30mm | Order number: OCZ/30 | Additional access port 30 mm |
| 60mm | Order number: OCZ/60 | Non-standard access port 60 mm |
| 100mm | Order number: OCZ/100 | Non-standard access port 100 mm |
| | Order number: KD | Access control |
| | Order number: CCTV | Camera inside the chamber |
| • | Order number: CCTV CAM | Additional camera |

