

# Data sheet

# Cooled Incubator (ST) ST 700 Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

#### Advantages of the SMART controller:

- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct downloadQuick change of program parameters
- Alarm Bar
- operating with gloves on



Smart - preview screen



tel. (+48) 32 453 91 70




TECHNICAL DATA			
air convection	forced		
chamber capacity [I]	625		
working capacity [I]	611		
controller	microprocessor PID		
display	4,3" full colour touch screen		
TEMPERATURE			
temperature range [°C]	+3+40 / +70 (opcja)		
temperature resolution every [°C]	0,1		
temperature fluctuation at 37°C [+/-°C]*	0,3		
temperature variation at 37°C [+/-°C]*	1,0		
over temperature protection	class 1.0 to DIN 12880 / class 3.3 (option)		
CHAMBER			
door type	solid / glass or double (option) /4/		
interior			
B Smart	aluminium		
C Smart	stainless steel to DIN 1.4016		
CS Smart	stainless steel to DIN 1.4016		
P Smart	DIN 1.4301		
PS Smart	DIN 1.4301		
housing			
B Smart	powder coated sheet		
C Smart	powder coated sheet		
CS Smart	stainless steel polished		
P Smart	powder coated sheet		

stainless steel polished

overall	dims	[mm]	/1/

PS Smart

overall dillis [illili] /1/	
width A	750
height B	1990
depth C	890
internal dims [mm]	
width D	540
width D'	600
height E	1510
depth F	680
height I	1380

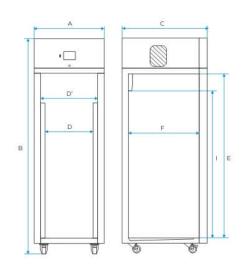

shelves (standard   max)	3   11
max shelf workload [kg] /2/	30
- reinforced shelf version (PW) [kg] /3/	100
max unit workload [kg]	150
weight [kg]	121





A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C

ELECTRICAL PARAMETERS		
voltage**	230V 50-60Hz	
nominal power [W]	650	
refrigerant	R290 / GWP=3	
warranty	24 months	
manufacturer	POL-EKO-APARATURA	

all the above technical data refer to standard units (without optional accessories)

\* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

K = +/- (T average max. - T average min.) / 2

- \*\* other power supplies on request
- 1 depth doesn't include 50 mm of power cable
- 2 on uniformly loaded surface
- 3 reinforced shelf
- 4 additional internal glass door

## **OPTIONS AND ACCESSORIES**

	Order number: */C	Internal glass door
	Order number: */A	External glass door
	Order number: */P	Wire shelf
INOX	Order number: */P INOX	Stainless steel wire shelf INOX
	Order number: */PP	Perforated shelf
Alla		



Order number: KUW GN\*/\*

Stainless steel cuvettes



Order number: ST/CHL/SWP ALU Aluminum drawer with powder coated slides

Reinforced shelf



Order number: ST/CHL/SWP INOX Stainless steel drawer with powder coated slides



Order number: ST/CHL/SWPN INOX Stainless steel drawer with stainless steel slides



Order number: \*/PW



A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C



Order number: QLK\* Castors



Order number: ST/70 Extended temperature range to 70°C



Order number: \*/FOT Photoperiodic system



Order number: GNZ Internal socket



Order numberLabDesk LabDesk software



Order number: BRT/\*/L or IQ/OQ/PQ Calibration and IQ, OQ, PQ qualification

DIN 3.3

Order number: \*/3.3 Over temperature protection 3.3 class according to DIN 12880



Order number: \*/RK Chart recorder



Order number: ZKM Magnetic door lock



Order number: BPP 12 Battery backup for display



Order number: PORT ALARM Dry alarm contact



Order number: USBK USB cable



Order number: LANK LAN cable



Order number: OCZ/20 Non-standard access port 20 mm



Order number: OCZ/30 Additional access port 30 mm



Order number: OCZ/60 Non-standard access port 60 mm



Order number: OCZ/100 Non-standard access port 100 mm





tel. (+48) 32 453 91 70

e-mail: info@pol-eko.com.pl internet: www.pol-eko.com.pl

### POL-EKO-APARATURA SP. J.

A. POLOK – KOWALSKA, S. KOWALSKI 44-300 WODZISŁAW ŚLĄSKI, UL. KOKOSZYCKA 172 C



Order number: KD

Access control

