

SNOL 7.2/1300 ashing furnace with solid ceramic chamber is designed for ashing and burn off processes up to 1300°C. Fan-assisted chimney permits to eliminate smoke from chamber during processing.

DESCRIPTION

- ✓ Solid ceramic chamber,
- ✓ Partially exposed heating elements on four sides around the chamber,
- ✓ Outside casing – metal sheet, powder painted grey,
- ✓ Fan-assisted chimney for smoke and humidity extraction,
- ✓ Door opens sideways,
- ✓ Door safety interlock switch,
- ✓ Control panel is placed in the underpart of the furnace,
- ✓ Non-programmable temperature controller – Omron E5CC,
- ✓ Ceramic bottom plate,
- ✓ Low power consumption,
- ✓ Good stability and uniformity,
- ✓ 1 year warranty.

Technical data	Units	Specifications
Useful volume	Liter	7.2
Rated power not more than	kW	3.5
Rated supply voltage	V	230
Rated frequency	Hz	50
Number of phases	-	1
Continuous operating temperature	°C	T+10-1300
Maximum temperature	°C	1300
Working chamber material	-	Ceramic
Maximum heating- up time (without charge),	Min.	240
Temperature stability in working chamber at rated temperature in thermal steady state without charge not more than	± °C	2
Temperature uniformity in working space at rated temperature in thermal steady state without charge not more than	± °C	10
Furnace working chamber dimensions:		
width	mm	195
depth	mm	295
height	mm	120
Furnace dimensions:		
width	mm	645
depth	mm	710
height	mm	705*
Mass (Netto)	kg	103

* – should increase by 20-85 mm, depending on the chimney

CONTROL:

- ✓ Temperature measurement by thermocouple type "S".
- ✓ PID electronic regulator, double digital display reference temperature and measured temperature.
- ✓ Products are equipped with high-precision digital microprocessor Omron or Eurotherm temperature controllers fitted with self-tuning and manual PID settings. The customer can select a basic or programmable temperature controller, which offers up to 32 programming segments (rate of temperature rise or decrease control, maintenance of present temperature, automatic shutdown). A wide range of devices allows selecting the most appropriate controller for your process.
- ✓ SSR control unit.

SUPPLIED DOCUMENTS:

- ✓ Furnace and temperature controller instructions,
- ✓ Electric diagram

OTHERS:

- ✓ CE marked

PACKING:

- ✓ Wooden box

OPTIONS:


- ✓ Eurotherm 3216 (non-programmable)
- ✓ Eurotherm 3208 (programmable)
- ✓ Omron E5CC-T (programmable)
- ✓ PC connection and SNOL software
- ✓ OTP (over temperature protection, non-adjustable)
- ✓ OTP (over temperature protection, adjustable Eurotherm 3216i)
- ✓ OTP (over temperature protection, adjustable Omron E5GC)
- ✓ Buzzer
- ✓ Timer (delayed furnace start)
- ✓ Ventilation hole ceramic tube (in the back side Ø17mm)

WARRANTY:

- ✓ One year limited warranty and later service for furnace
- ✓ Possibility to extend warranty for an additional 1 year

COUNTRY OF ORIGIN:

- ✓ Lithuania (EU country)

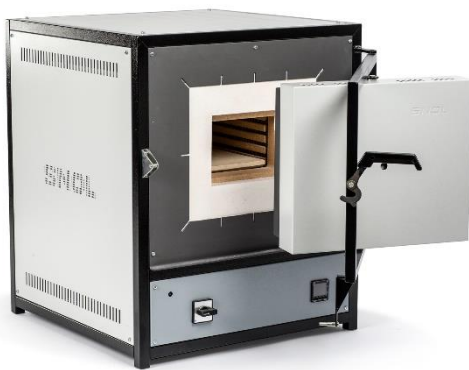

Administrative address
Umega Group, AB
Metalo str.5, 28216 Utena
Lithuania

Factory address
Umega Group, AB, SnolTherm unit
Plento str.3, 28104 Utena
Lithuania

Tel. +370 389 54586
sales@snoltherm.com
www.snol.com

VAT code: LT263347219
Company registration No. 126334727