


PrO-Cyt

By Centurion Scientific Ltd



LabTeam Scandinavia AB
Tele Sverige: 042-300 91 30
Tele Danmark: 36 98 02 30

www.labteamet.com
info@labteamet.com
order@labteamet.com

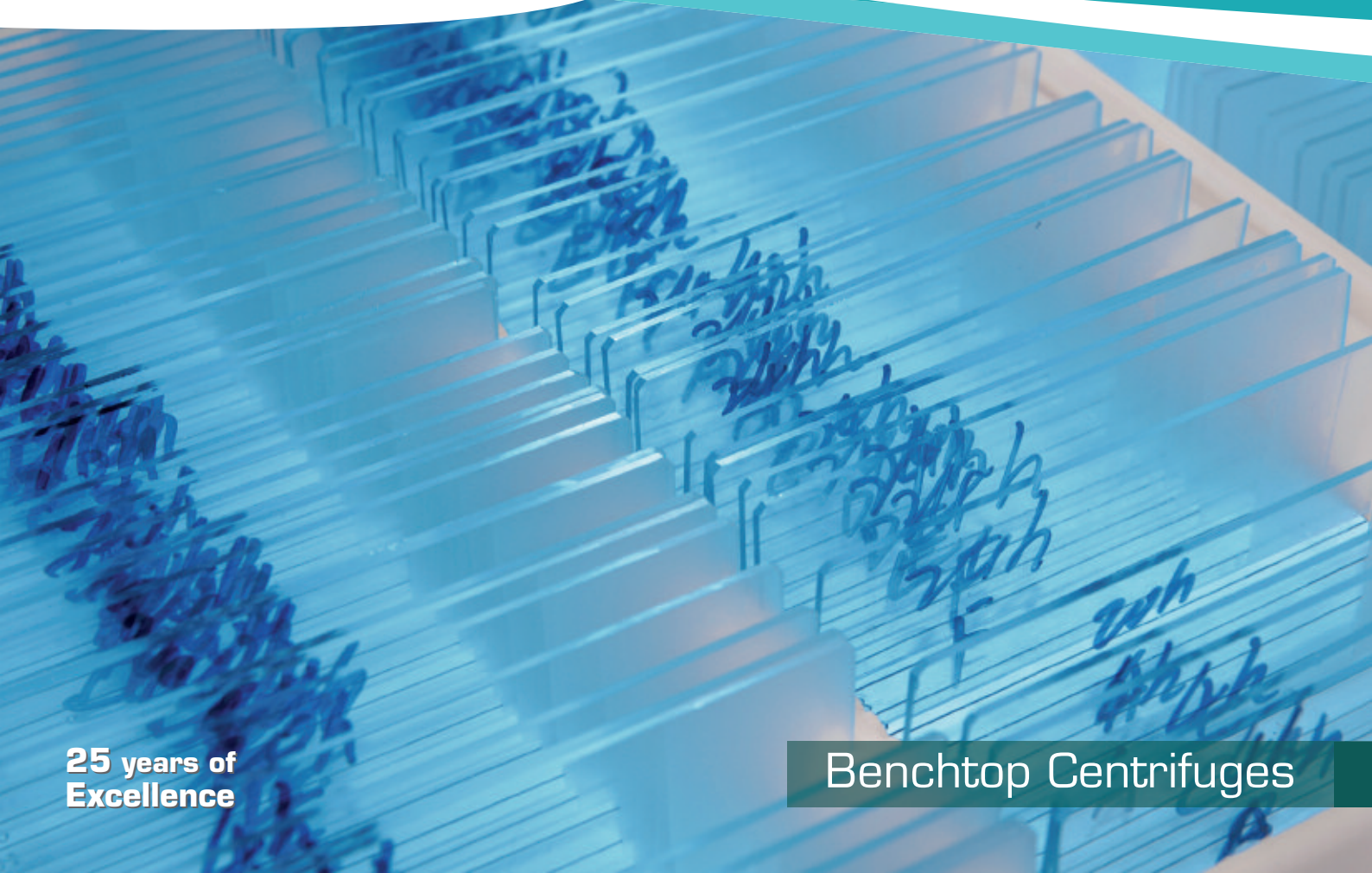

**25 years of
Excellence**

Benchtop Centrifuges



CE marked and ISO 9001
regulated and now IVD certified too

The PrO-Cyt Range

Our PrO-Cyt range are compact and practical centrifuges for small , medium and large sample volumes. Unlike others we offer you a choice of two Centrifuges and each can use one of three rotors , from 4 , 8 and 12 places , thus meeting your needs exactly.

They are both cost effective, yet still benefiting from the superb build quality and engineering of our research series.

Centurion centrifuges, with their applications-oriented accessories, are proven favourites in clinical, cytology , histology and small medical laboratories. Our compact, ergonomic design: smaller dimensions, ergonomic height and curved features ensure a comfortable and efficient workspace. Our Centrifuges deliver optimum performance in compact units, designed for ease of use. Quiet (below 50Db, rotor dependant) with low heat operation noise levels have been reduced while maintaining sample integrity with good airflow .

Performance and safety are top of our agenda, tight operational control is maintained by the microprocessor. An imbalance sensor is incorporated into the control loop to bring to a controlled stop of the centrifuge in the case of a tube rupture or incorrect loading.

Centurion centrifuges have been an industry standard for quality and reliability for 25 years, our centrifuges can be relied upon to deliver safe, reproducible and trouble-free results for many years.

Warranty 3 years for 2014 models

Centurion brand centrifugation products are part of a dedicated range of products and services designed to provide you with quality equipment and support.

Choosing The Right Centrifuge

A centrifuge is a critical piece of equipment. Therefore, when the time comes to add new centrifuges or replace old models, it is important to make an informed decision. Bringing together technological and design excellence, Centurion have paved the way for a new generation of centrifuges.

User features to PrO-Cyt range	Advantages
Large LCD display with touch control	Clarity & ease of use
10 acceleration rates 10 deceleration rates	Sample assistance
108 memory (12 pages of 9)	Multi department & users
Timer 0-99 minutes & Hold in 10 second increments	Total flexibility
Pulse short run	Fast pelleting
Run in Speed (Rpm) or Rcf (G) in 1 Rpm increments	Accuracy
Sound <55Db	Quiet
Orientation acceleration rate	Prevents initial sample side deposits

Standard to PrO-Cyt range	Advantages
Extra thick Stainless steel bowl	Easy clean & rust free
Port to lid	Tachometer speed certification
Alloy & steel frame with composite's	Strong light & quiet
World leading Industrial grade Inverter	Reliability & strength
Best Quality European Brushless motor	Quiet, cool & reliable
High technology airflow (ambient models)	Cooler running
All centrifuges have user accessible service sections to access all safety parameter	

Safety features to PrO-Cyt range	Advantages
Multi point lid locking	Lid safety
Emergency lid release	Power failure release
Lid gas strut	Lid dropping safety
Lid lock detection	To run Lid must be shut
Imbalance detection	Easy "eye" volume loading
Overspeed Sensor	Safe detection of speed
Set inverter values	Electronic safety of speed
Barrier ring	Extra metal protection of chamber
Motor overheat sensor	Safe motor protection

Centurion Scientific Ltd centrifuges comply to all relevant EU standards for Quality and medical Devices IEC 61010 and CE conformity test marks Emission, immunity to EN/IEC 61326-1, Class B

Index

Page 4 PrO-CytLCD Centrifuges

Page 5 Rotors available. 4 place , 8 place and 12 place

Page 6 Sample holders available, single, double or single large.

Page 7 Package supplied and other available items

PrO-CytLCD



Display indicative only.

- Speed** 200-2000 Rpm
(1 Rpm steps)
- Rcf max** 550 G
- Timer** 0-9999 Mins & hold
(1 sec steps)
- Dimensions** ... HWD 275 x 325 x 470mm
- Weight** 21.8 Kg (without rotor)
- Power** 130 watts
- Memory** 108 programs
(12 pages x 9)
- Accel rates** ... 10 programs
- Decel rates** ... 10 programs

Choosing The Right Centrifuge

Pro-CytLCD and its Cytology rotors are 100% compatible with Shandon Cytospin®, this allows you to continue using exactly the same protocols and consumables as on a Cytospin®.

If you wish to buy your cyto- consumables directly from us, please contact our product specialist that will provide you information and prices.

Available Cytology Rotors

Using centrifugal force to separate and deposit a monolayer of cells onto slides whilst maintaining cellular integrity within a clearly defined area FROM ANY FLUID MATRIX. Cyto centrifugation also constructively flattens cells for excellent nuclear presentation.

Applications include Cytology, Histology, Haematology, Oncology, Immunochemistry, Virology, Serology and Microbiology.

Offering samples from 0.1ml to 6ml and an extensive range of accessories we have your Cytology needs covered.

Cytology rotor includes Rotor, sealed lid, Clips, 24 x holders with filter card & 100 slides

Rotors	4420
Rotor type	4 x 0.2ml
Holders	Single or double
Minimum speed	200 Rpm
Maximum speed	2,000 Rpm
Maximum Rcf (G)	550
Radius max cms	12
Sample tube angle (o)	0
Acceleration time (secs)	20
Deceleration time (secs)	20
Autoclaveable (frequency)	121°C (10)



Rotors	4430
Rotor type	8 x 0.2ml
Holders	Single or double
Minimum speed	200 Rpm
Maximum speed	2,000 Rpm
Maximum Rcf (G)	550
Radius max cms	12
Sample tube angle (o)	0
Acceleration time (secs)	20
Deceleration time (secs)	20
Autoclaveable (frequency)	121°C (10)



Rotors	4460
Rotor type	12 x 0.2ml
Holders	Single or double
Minimum speed	200 Rpm
Maximum speed	2,000 Rpm
Maximum Rcf (G)	550
Radius max cms	12
Sample tube angle (o)	0
Acceleration time (secs)	20
Deceleration time (secs)	20
Autoclaveable (frequency)	121°C (10)



Accessories

4444

Single sample holder
with card (up to 1ml)



4446

Double sample holder
with card (up to 1ml)



4600

Single sample holder
with card (up to 6ml)



4462

Stainless steel clips



Accessories included

All Pro-Cyt centrifuges are complete with a starter pack this includes.

100 microscope slides and a pack of 24 double sample holders including fixed blotters

Rotors are complete with lid , swing plate and stainless steel holders to match.



We trust that the details contained within these pages allows you to make a informative decision on our Pro-Cyt range.

If however you need more information or assistance, please do not hesitate to call