

himac centrifuges

ISO 9001 BUREAU VERITAS Certification

Centrifuge & Processing Equipment Division of Koki Holdings Co.,Ltd. is certificated with ISO 9001, Quality management system.

CAUTION: • For safe and proper use of your machine, carefully read and follow instruction manual.

• All specifications are subject to change without advanced notice.

- Actual color may vary from the color the photos on this catalogue, due to printing condition.
- Due to safety reason, installation environments, operating environments and conditions may be restricted.
 Unless specially mentioned, products and/or operation panel of the photos are standard specifications.
- CT18R centrifuge is not a medical device.
- For further information. contact your nearest sales representative.

himac

Koki Holdings Co., Ltd. Centrifuge & Processing Equipment Division Shinagawa Intercity Tower A, 2-15-1 Konan, Minato-ku, Tokyo 108-6018 Japan Tel : +81-3-6738-0863 Fax : +81-3-6738-0861 https://www.himac-science.com




Multipurpose Centrifuge



Operation Panel has 2 status indicators!

[Status 1]

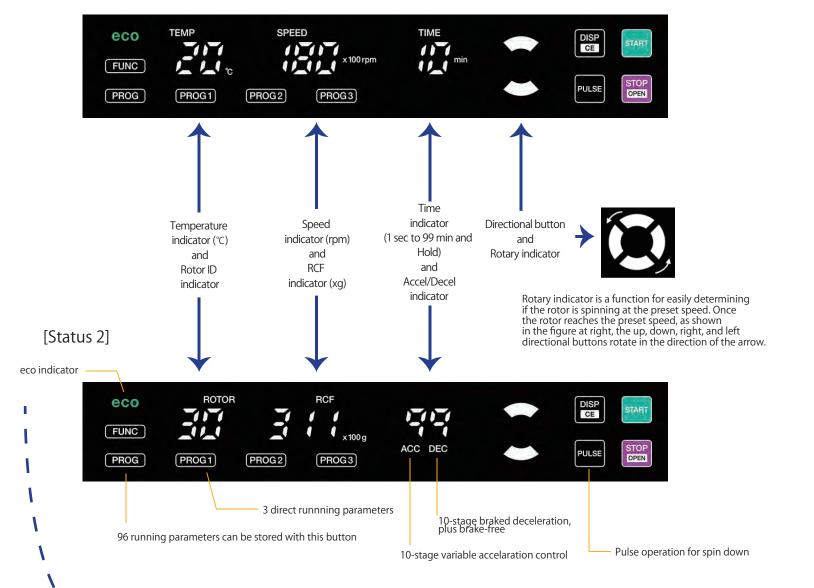

Actual Run Time (ART Function)


Timer starts counting after the set speed is attained. So acceleration time is excluded from the centrifugation time. You can obtaion higher reproducible centrifugation result with this function. It easier to obtain more consistent separation results because there is less variation in acceleration times due to different rotors and sample volumes.

Total Run Hours Indication

Total run hours^{*1} can be displayed on the operation panel. It is useful as a guideline of the periodical maintenance. ^{*1} Including acceleration time regardless a setting of the Actual Run Timer

Optional Remote monitoring Sytem himac View

Remote monitoring application "himac View" is available for CT18R, CR30NX and CP-NX series. "himac View" monitors, operates of start or stops a centrifuge with iPhone™ or Android ™ Phone via Wi-Fi®. "himac View" monitors up tp 16 units of CT18R, CR30NX or CP-NX series in the LAN access network. * Optional LAN board in CT18R is required to use "himac View". For more details contact the nearest sales representative.

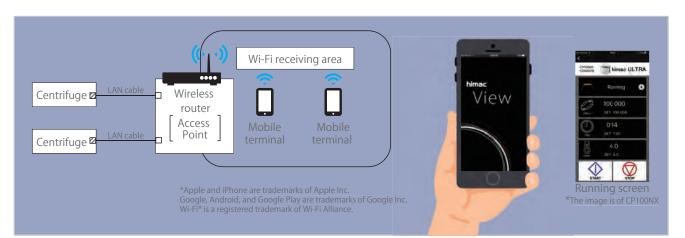

Self-Locking Rotor System

Just place a rotor on the drive shaft. The rotor is automatically locked by centrifugal force once it starts spinning. When the rotor comes to complete stop, the lock is automatically released.



Power Saving Function

When "eco setting" is enable, if the door is not opened or closed and no operations are performed for 10 minutes, the cooling device is automatically stopped.

Max. 25% stand-by electricity is reduced,

- comparing the stand-by electricity without eco setting. (under in-house
- evaluation) When a user forgets to turn off the main power switch

for a long time (such as overnight), stand -by electricity consumption is

reduced.

Set "eco" function and contribute energy saving!



Secure Door-Lock system: Press the door down softly then the door is locked automatically (Both right side & left side)

Rotors, Buckets, Adapters, Tubes & Bottles

Fixed Angle Rotors Auto clave Auto clave T18A41 Auto clave T15A46 T15A39 0 The rotor marked with this icon can be sterilized by autoclave at 121°C for 20 minutes. Auto clave 36 x 2.0/1.5ml, 18 x 0.2ml 18,000rpm (31,100xg/23,200xg) 24 x 2.0/1.5ml 15,000rpm (21,500xg) 12 x 5ml 15,000rpm (22,000xg) * Rotor cover is sold separately. T15A41 T12A42 T12A41 Auto T9A31 Auto 4 x 15TC/50TC 6 x 50TC 6 x 80ml 8 x 50TC 15,000rpm (22,260xg) 12,000rpm (18,000xg) 9,000rpm (10,600xg) 12,000rpm (17,100xg)

Rotor	Rotor	Max.	Max.	Nominal	Tube Adapter		Tube / Bottle				
Model (Part No.)	ID	speed (rpm)	RCF (xg)	Capacity (pcs x ml)	P/N	Name	P/N	Name	Tube Size (\$\phi x L) mm	Se ll ing Unit	
T18A41 (91570000)		18,000	Out: 31,100 In: 26,400	36 x 2/1.5	_		_	• 2.0ml/1.5ml Microtube	10.8 x 39	_	
	30	18,000	Out: 27,100 In:22,400	36 x 0.7		-	_	• 0.7ml Microtube	10.8 x 25	—	
		15,100	Out: 19,900 In: 16,600	36 x 1/0.5	S413789A	1GA2 Adapter (2pcs/set)	489241A 489240A	1.0ml Glass Tube (U-bottom) 0.5ml Glass Tube	9.0 x 35 9.0 x 35	100pcs/set 100pcs/set	
		15,100	Out: 20,800 In: 17,500	36 x 0.5	S414002A	0.5 Adapter (B) (4pcs/set)	489239A	0.5ml Glass Tube (U-bottom)	6.0 x 47	100pcs/set	
		18,000	Out: 28,500 In: 23,800	36 x 0.5	S414004A	0.75 Adapter (B) (4pcs/set)		• 0.5ml Microtube	7.8 x 30	_	
		18,000	Out: 30,100 In: 25,400	36 x 0.4/0.25	S414000A	0.4 Adapter (B) (4pcs/set)		 0.4ml Microtube 0.25ml Microtube 	5.7 x 46 5.7 x 30	-	
		18,000	23,200	18 x 0.2				0.25ml Microtube 0.2ml PCR Tube	6.4 x 24	_	
		15,000	21,500	24 x 2/1.5	_		-	2.0ml/1.5ml Microtube	10.8 x 39	—	
	49			24 x 0.7		_	_	 0.7ml Microtube 	10.8 x 25	_	
				24 x 1/0.5	S413789A 1GA2 A	1GA2 Adapter (2pcs/set)	489241A	1.0ml Glass Tube (U-bottom)	9.0 x 35	100pcs/set	
T15A39						IGAZ AUZPIEL (ZPCS/SEL)	489240A	0.5ml Glass Tube	9.0 x 35	100pcs/set	
(91589900)				24 x 0.5	S414002A	0.5 Adapter (B) (4pcs/set)	489239A	0.5ml Glass Tube (U-bottom)	6.0 x 47	100pcs/set	
				24 x 0.5/0.4	S414004A	0.75 Adapter (B) (4pcs/set)	_	• 0.5ml Microtube	7.8 x 30	_	
				24 x 0.4/0.25	S414000A 0.4 Adapter (B) (4pcs/set)		0.4ml Microtube	5.7 x 46	—		
						—	0.25ml Microtube	5.7 x 30	_		
T15A46 (91583600)	39	15,000	22,000	12 x 5		_	Eppend	orf Tubes [®] 5.0mL with screw cap)		
		15,000	0 22,260	······································		S307904A	himac 50TC tube	30 x 115	300pcs/box		
T 45A 44	42	15,000	22,200	4 x 15 (TC)		_		himac 15TC tube	17x120	100pcs/box	
T15A41 (91591300)		15,000	14,200	4 x 5	 Universal adapter by Eppendof AG 		Eppendorf Tubes [®] 5.0mL with screw cap				
		11,000	11,970	4 x 50 (TC) 4 x 15 (TC) or		_		● 50ml TC Tube (V-bottom) ● 15ml TC Tube (V-bottom)	30 x 115 17 x 120		
T12A41 (91570100)	20	12,000	18,000	6 x 50 (TC)	-		S307904A	himac 50TC tube	30 x 115	300pcs/box	
				6 x 50 (TC)	_		—	 50ml TC Tube (V-bottom) 	30 x 115	—	
		10,000	12,200	6 x 15 (TC)	S411576A 15TA50T adapter (C) (2pcs/set)			himac 15TC tube		100pcs/box	
					S411576A 1	5TA50T adapter (C) (2pcs/set)	_	● 15ml TC Tube (V-bottom) ● 50ml TC Tube (V-bottom)	17x120		
T9A31 (91587400)	24	9,000	10,600	8 x 50 (TC) 8 x 15 (TC)			_	● 15ml TC Tube (V-bottom)	30 x 115 17 x 120	_	
			<u>10,300</u> 9,700	8 x 15 (1C) 8 x 13		3A50T adapter (B) (2pcs/set)	_	 13ml Tube (U-bottom) 	17x120 17x101	-	
T12A42 (91570200)	21	12,000	17,100	6 x 80			332094A	80PA Bottle	38x106	10pcs/set	
				6 x 70		_	S308433A	70PC Bottle (B)	38x106	6pcs/set	

Swinging Bucket Rotors



Rotor	Rotor	Max.	Max.	Nominal	Tube Adapter		Tube / Bottle				
Model ID	ID	speed (rpm)		Capacity (pcs x ml)	Bucket	P/N	Name	P/N	Name	Tube Size (Ф x L) mm	Selling Uni
		4,200	2,620	32 x 5		S307332A	5x8S Adapter (4pcs/set)	• 4-5ml Bl	ood Collection Tube	12.4-13.2 x 75-78	—
		4,200	3,160	2210		S307333A	10x8S Adapter	7ml Bloo	d Collection Tube	12.4-13.2 x 100	—
		3,600	2,320	32 x 10		S307333A (4pcs/sel)	464977A	10 Glass Tube	12 x 105	50pcs/box	
		4,200	3,160	- 32 x 15	Standard accessory	S307334A	15x8S Adapter (4pcs/set)	10ml Blo	od Collection Tube	15.4-16 x 100	_
		3,600	2,320	JZ X 15				464976A	15 Glass Tube	16.5 x 105	50pcs/box
T4SS32	61	4,200	3,200	24 x 15		S307335A	15TCx6S Adapter (4pcs/set)	• 15ml TC	Tube (V-bottom)	17 x 120	_
(91570300)		4,200	3,160	8 x 50		S307336A	50TCx2S Adapter (4pcs/set)	• 50ml TC	Tube (V-bottom)	30 x 115	_
		3,500	2,120	8 x 50		S308844A	50x2S Adapter (4pcs/set)	473866A	50 Glass Tube	35 x 100	10pcs/box
		3,000	1,640	40 x 15		S308546A	15x10S Adapter	464976A	15 Glass Tube	16.5 x 105	50pcs/box
		4,200	3,210	40 x 10		0000010/1	(4pcs/set)	10ml Blo	od Collection Tube	15.4-16 x 110	—
		4,200	3,090	40 x 10		S308695A	10x10SB Adapter (4pcs/set)	• 5-10ml B	Blood Collection Tube	12.4-15.6 x 75-100	
P/N : S410523C Round Bucket (Option) (4pcs/set) Adapters and a cap (Option) for T5SS31 can be used.						1					
	50	4.800 4.170	4,170	28 x 5	-	S409812A	5x7S Adapter (4pcs/set)	• 4-5ml Blood Collection Tube 12.4-13.2 x 75-78		-	
		4,000	4,170	28 x 7		S409811A	10x7S Adapter (4pcs/set)	7ml Bloo	d Collection Tube	12.4-13.2 x 100	—
		3,000	1,630	28 x 10			(4pcs/set)	464977A	10 Glass Tube	12 x 105	50pcs/box
		4,800	4,170	28 x 10		S409810A 1	15x7S Adapter (4pcs/set)	10ml Blo	od Collection Tube	15.4-16 x 100	—
		3,000	1,630	28 x 15				464976A	15 Glass Tube	16.5 x 105	50pcs/box
			3,990	4 x 50		S409813A	50x1S Adapter (4pcs/set)	473866A	50 Glass Tube	35 x 100	10pcs/box
T5SS31		4,800 4,170	4,170	12 x 15(TC)		S409976A	15TCx3S Adapter (4pcs/set)	Optional bu	Tube (V-bottom) icket cap can be used.	17 x 120	—
(91587700)				16 x 15(TC)		S409814A	15TCx4S Adapter (4pcs/set)	 15ml TC Optional but 	Tube (V-bottom) ucket cap cannot be used.	17 x 120	—
				4 x 50(TC)		S409815A	50TCx1S Adapter (4pcs/set)	• 50ml TC	Tube (V-bottom)	30 x 115	—
				4 x 160		S410251A	160x1S Adapter (4pcs/set)	S306847A	160PP Bottle**	52 x 120	10pcs/box
						011020111	(4pcs/set)	S306848A	160PC Bottle**	52 x 120	10pcs/box
				4 x 250		_		250PF	PCO Bottle	60.9 x 133	Nalgene [™] P/N 3141-02
			4 X 2 J 0				• 250PC	Bottle	61.5 x 135	Nalgene [™] P/N 3140-02	
		P/N : S410	252A Shield C	ap (Option) (4pc	s/set) The cap cannot	be used for	15TCx4S Adapter	and 250ml	Bottle.		
T5S33 (91591700)	56	4,700 3		4 x 96DW		_		96-hole Deep Well Plate (128x86x43H)			
			3,100	12 x 96MP	Standard accessory			96-hole Micro Plate (128x86x14H)			
			2 x F					 Deep Well Plate with Filter (128x86x74H) 			

Remarks

Items with ● are marketed products. Purchase directly from your local supplier.
 Actual volume of the tube or bottle may be less than the nominal volume due to design of the rotor or shape of the tube/bottle.
 SOmL TC tubes and 15mL TC tubes cannot be centrifuged together in T15A41 rotor.
 ** Cap of 160PP bottle and 160PC bottle are sold separately. Please order P/N: S410008A 160 Cap (2pcs/set) with the bottle.
 When "● 50ml Bottle (U-bottom)" is used, please refer to data sheet of the bottle and set RCF at Max.RCF or lower.
 Eppendorf Tubes® is a registered mark of Eppendorf AG.
 Nalgene™ is a trademark of Thermo Fisher Scientific Inc.

High performance himac TC Tubes!

himac original 50 ml & 15ml conical tubes (centrifugation at high RCF)

- (1) 50ml tubes can be used up to max. 50,000 ×g *[actual capacity, 45 ml/tube]
- (2) 15ml tubes can be used up to max. 36,100 ×g **
- (3) V-shaped bottom collects sediments more easily at the bottom of the tube compared to rounded-bottom tubes.
- (4) Having an air-tight screw cap. Improved safety!
- (5) Radiation Sterilized
- (6) Volume scales and space for labels

* with R19A2 rotor for model CR30NX ** with R16A2 rotor for model CR30NX



S311184A: himac 15TC TUBE (100pcs/set) S307904A: himac 50TC TUBE (300pcs/set)

Centrifuge Model	Rotor Model	TC Tube Type	Max. Speed (rpm)	Max. RCF (xg)	No. of tubes X Actual capacity (ml)	Required Adapter
CT18R	T15A41	himac 50TC Tube	15,000	22,260	4 X 45	_
		himac 15TC Tube	15,000	22,260	4 X 15	_
CR22N	R18A	himac 50TC Tube	18,000	42,200	8 X 45	-
Max		himac 15TC Tube	10,000	12,700	8 X 15	S411576A (2pcs/set)
	R15A	himac 50TC Tube	15,000	32,300	10 X 42	_
	E	himac 15TC Tube	15,000	32,500	10 X 15	-
CR30NX	R19A2	himac 50TC Tube	19,000	50,000	8 X 45	_
		himac 15TC Tube	9,650	12,610	8 X 15	S411576A (2pcs/set)
	R16A2	himac 50TC Tube	16,000	36,100	10 X 42	_
		himac 15TC Tube	16,000	36,100	10 X 15	_



Maximum speed	18,000 rpm (1
Maximum RCF	31,100×g (T
Maximum capacity	1,000 ml (4 X
Set speed range	300 rpm to 18
Temperature range adjustment	-20°C to +40°
Timer	1 to 99 min. (in 1 to 59 sec. wit
Acceleration/Deceleration Mode	Accel: 10 mod
Memory-based programmed operation	99 programs,
Drive Unit	Induction mot
Refrigerator	300W Herme
Rotor stop signal	Select from 5
Interface connection	LAN X 1 (Optio
Safety functions	Door interlock, imbalance dete
Noise level	Less than 55dB Less than 59dB
Compliance with standard	CE, cCSAus
Heat radiation	650W
Power requirements	AC220-240V ± AC115-127V ±
Dimensions	634 mm (W) 3
Weight	80 Kg
Installation Environment	Ambient tempe Ambient tembe
	* for the United **Ordering num







Top View

Note: R18A, R15A, R19A2 and R16A2 have a rotor cover as a standard accessory.

For more details of CR22N and CR30NX, visit our website: https://www.himac-science.com

(T18A41) T18A41) (250 ml by T5SS31) 18,000 rpm °℃ in increments of 1 min.) with HOLD function or ith HOLD function odes, Decel: 10 modes and brake-free , including 3 direct select buttons otor (Inverter) etic compressor (Refrigerant: HFC 134a) 5 melodies, beep and mute ion) k, dual overspeed detector, Rotor ID detector, tector, over-temperature detector IB by fixed angle rotor B by swinging bucket rotor ±10% 50Hz 8A (ordering number: 915617C1 and 915617C2**) ±10% 60Hz 15A (ordering number: 915617C3)

X 550 mm (D) X 340 mm (H) (Height to the table: 291 mm)

erature for operation:2°C to 40 °C erature for assuring performance: 15℃ to 25℃

d States and Canada: AC115-127V ±10% 60Hz 15A (ordering number: 915617C3), only *Ordering number 915617C2 is applicable to EU region