





LF 1400

HIGH LEVEL OF SECURITY FOR YOUR SAMPLES, DUE TO THE GREAT COOLING PERFORMANCE WITH AIR FORCED COOLING.

Automatic defrost saves you time to focus on what's most important and stores the samples safely.



AVAILABLE IN

ATEX





AIR FORCED FILT COOLING CONS

FILTERLESS CONSTRUCTION

STAINLESS STEEL & SPARK FREE

The unit is available with stainless steel AISI 304 Exterior, so that the interior and exterior is made entirely of stainless steel, providing a stronger more resistant unit.

This unit can be made in an ATEX version – approved according to EN 60079-15 category 3, zone 2. CE EX 11-/3 G

FLEXIBLE STORAGE SOLUTIONS


Arctiko's LF range is a highly customizable allowing you to choose between a wide range of universal shelves and drawers.


This allows the user to organize their unit exactly as desired. Our universal shelf and drawer guides allow you to mix shelves and drawers for flexible storage.


GREAT COOLING PERFORMANCE

Air forced cooling provides superior temperature uniformity and stability. Forcing the airflow around the shelves and the stored items, ensures that all parts of the unit is equally cold and avoids any cold or hot spots. Distributed by a powerful fan and open air flow areas along the sides, the air is always distributed evenly throughout the fully loaded unit.

MAIN DIMENSION FOR LF 1400

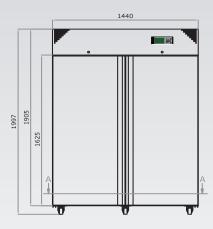

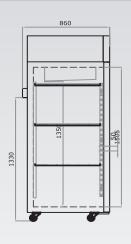

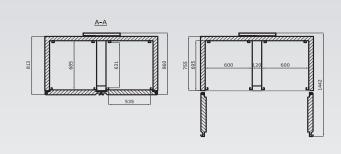

MAIN SPECIFICATIONS

Item code	DAI 0275 / DAI 0276
Temperature range (°C)	-10 / -30
Max. ambient temperature (°C)	25
Capacity (L)	1361
Exterior dimensions WxDxH (mm)	1440x860x1997
Interior dimensions WxDxH (mm)	1320x685x1505
Weight (KG)	195
Insulation (mm)	60
Cooling technology	air forced
Controller model / probe type	G214 / PT1000
DOWED	

POWER 230' Power supply (V) Frequency (Hz) 50 Power consumption (kWh/24H) 18.1 AMP (A) 4.9 Power (Watt) 964 Horse power (HP) 1 1/4 ALARM FUNCTIONS Alarm display as text - not codes yes Visual / acoustic alarm yes / yes Power failure alarm yes Adjustable high / low temperature alarm yes / yes Open door alarm yes Probe failure alarm yes Low battery alarm yes Compressor failure alarm no Battery backup for alarms (approx. hrs.) 72 DATA LOGGING & EXTERNAL CONNECTION Voltage free contact for remote alarm (e.g. GSM alarm module) yes Electronic data logger ves RS 485 port yes USB port for software update & data download yes Temperature chart recorder optional DISPLAY FEATURES Display type 4,3" with touch buttons Number of optional reference probes 2 Password protection for turning unit on / off yes Temperature graph on display yes OTHER FEATURES Automatic defrost yes Antifreeze function no Vacuum release port no Access port for external probe etc. (int. diameter mm) 12 Light no Lock with key yes Castors (total pcs) / castors with brake (pcs) 5/3 Fixed feet (pcs) / levelling feet (pcs) no Sub lids / inner doors (pcs) no

Sub lids / inner doors (pcs)noReversible door(s)yesFixed shelves (pcs) / adjustable shelves (pcs)0 / 6Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)46 / 22Max. load per shelf (kg)30Adjustable drawers (pcs)optionalMax. no. of adjustable drawers (pcs)16Max. load per adjustable drawer (kg)30REFRIGERATION & CABINET

Number of compressors 1 Refrigerant R290 Interior cabinet material stainless steel Exterior cabinet material painted steel / stainless steel Colour (exterior cabinet) white / stainless steel Foam type polyurethane VIP vacuum insulated panels no SHIPPING Shipping weight (KG) 278 Shipping dimensions WxDxH (cm) 156x101x220 Shipping volume (m³) 3,47

*1110V/60Hz available as option

Arctiko A/S

Drawers (stainless steel)

Modbus adapter Voltage stabilizer Voltage converters LED light

GSM alarm module

OPTIONAL ACCESSORIES

Temperature chart recorder Shelves (perforated, wire plastic coated, wire stainless steel)



