





# PR 1400

OFFERS A GREAT OVERVIEW OF THE STORED ITEMS DUE TO THE DOUBLE PANE GLASS DOOR. Designed to meet the demands of the

pharmaceutical & biotechnology industries providing reliability, accuracy, and sample security.







AVAILABLE IN

COOLING

FILTERLESS

### STAINLESS STEEL & SPARK FREE

The unit is available with stainless steel AISI 304 Exterior, so that the interior and exterior is made entirely of stainless steel, providing a stronger more resistant unit.

This unit can be made in an ATEX version – approved according to EN 60079-15 category 3, zone 2. CE EX 11-/3 G

### FLEXIBLE STORAGE SOLUTIONS

Arctiko's PR range is a highly customizable allowing you to choose between a wide range of universal shelves and drawers.


This allows the user to organize their unit exactly as desired. Our universal shelf and drawer guides allow you to mix shelves and drawers for flexible storage.


## **GREAT COOLING PERFORMANCE**

Air forced cooling provides superior temperature uniformity and stability even with frequently door openings. Forcing the airflow around the shelves and the stored items, ensures that all parts of the unit is equally cold and avoids any cold or hot spots. Distributed by a powerful fan and open air flow areas along the sides, the air is always distributed evenly throughout the fully loaded unit.


### **MAIN DIMENSION FOR PR 1400**



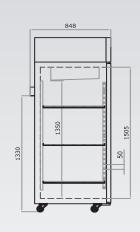

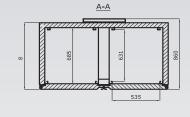



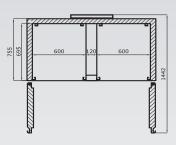

#### MAIN SPECIFICATIONS

Item code	DAI 0870 / DAI 0871
Temperature range (°C)	+1 / +10
Max. ambient temperature (°C)	32
Capacity (L)	1381
Exterior dimensions WxDxH (mm)	1440x860x1997
Interior dimensions WxDxH (mm)	1320x695x1505
Weight (KG)	209
Insulation (mm)	60
Cooling technology	air forced
Controller model / probe type	G214 / PT1000

Power supply (V)	230*
Frequency (Hz)	50/60
Power consumption (kWh/24H)	8,1
AMP (A)	3,8
Power (Watt)	618
Horse power (HP)	1/3

#### ALARM FUNCTIONS

Alarm display as text - not codes	yes
Visual / acoustic alarm	yes / yes
Power failure alarm	yes
Adjustable high / low temperature alarm	yes / yes
Open door alarm	yes
Probe failure alarm	yes
Low battery alarm	yes
Compressor failure alarm	no
Battery backup for alarms (approx. brs.)	72

#### DATA LOGGING & EXTERNAL CONNECTION

Voltage free contact for remote alarm (e.g. GSM alarm module)	yes
Electronic data logger	yes
RS 485 port	yes
USB port for software update & data download	yes
Temperature chart recorder	optional

#### DISPLAY FEATURES

Display type	4,3" with touch buttons
Number of optional reference probes	2
Password protection for turning unit on / off	yes
Temperature graph on display	yes

#### OTHER FEATURES

Automatic defrost	yes
Antifreeze function	no
Vacuum release port	no
Access port for external probe etc. (int. diameter mm)	12
Light	yes
Lock with key	yes
Castors (total pcs) / castors with brake (pcs)	5/3
Fixed feet (pcs) / levelling feet (pcs)	no
Sub lids / inner doors (pcs)	no
Reversible door(s)	yes
Fixed shelves (pcs) / adjustable shelves (pcs)	0/6
Max. no. of shelves (pcs) / recommended max. no. of shelves (pcs)	46 / 22
Max. load per shelf (kg)	30
Adjustable drawers (pcs)	optional
Max. no. of adjustable drawers (pcs)	16
Max. load per adjustable drawer (kg)	30

### REFRIGERATION & CABINET

Number of compressors	1
Refrigerant	R290
Interior cabinet material	stainless steel
Exterior cabinet material	painted steel / stainless steel
Colour (exterior cabinet)	white / stainless steel
Foam type	polyurethane
VIP vacuum insulated panels	no
SHIPPING	

Shipping weight (KG)	292
Shipping dimensions WxDxH (cm)	156x101x220
Shipping volume (m³)	3.47

<sup>\*110</sup>V/60Hz available as option

## OPTIONAL ACCESSORIES

GSM alarm module Temperature chart recorder Shelves (perforated, wire) Drawers (stainless steel) Modbus adapter Voltage stabilizer Voltage converters





